


TAG 3.0

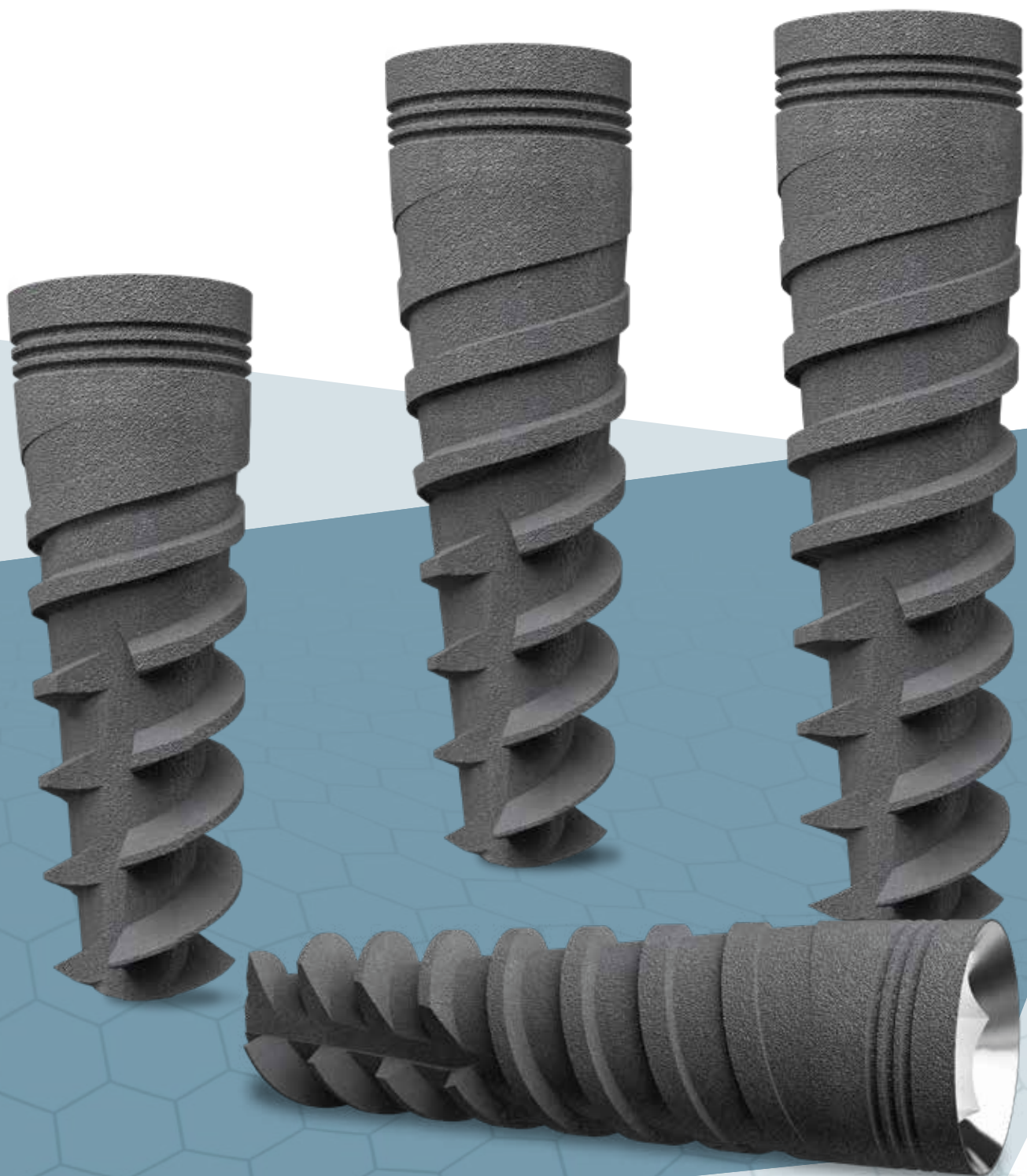

- Self-tapping conical dental implant
- Internal hexagon connection

GRADE 5 TITANIUM

TAG 3.0 comes from the need to solve atrophic and thin frontal crests clinical cases and as a valid solution in lateral agenesis.

Implant design is characterized by condensing conical body with osteotome effect and neck designed for reduce trauma to the crestal area.

The 1,6 mm connection screw with deep engagement ensures a precision fit connection between the prosthetic parts and implant.


Surface Treatment "Sandblasting and acid-etching" helps to obtain a micro-roughness which enhances the speed of the osseointegration process.

Excellent decontamination is performed by using an Argon Plasma reactor in a cleanroom.

	10 mm	11.5 mm	13 mm
			
Tag 3.0 codes	TAG MF FIXTURE 3.0 H10 cod. TAGMF029	TAG MF FIXTURE 3.0 H11.5 cod. TAGMF030	TAG MF FIXTURE 3.0 H13 cod. TAGMF031

Ø 3.0 Fixture

Ø 3.5 Fixture

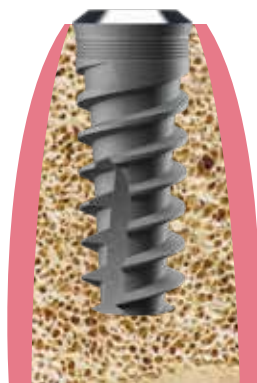

Thin crests



Narrow spaces

TAG 3.0 Prosthetic components

	Healing cap		Transfer open tray		Analog
	H.2 - cod. TAGVG32 H.4 - cod. TAGVG34		cod. TAGTRA3		cod. TAGANA3
	Castable cylinder N/R and Rot.		Straight abutment		Titanium temporary abutment N/R and Rot.
	cod. TAGCC3 cod. TAGCC3R		H.0 - cod. TAGMD3 H.2 - cod. TAGMD32		cod. TAGMP3 cod. TAGMP3R
	Angled abutment 15°		Angled abutment 25°		Abutment screw M1,6
	cod. TAGMA315		cod. TAGMA325		cod. TAGVM3
	Ball abutment		Open low cap with o-ring		Teflon cap cod. TAGCAPT
	H.1 - cod. TAGASF31 H.2 - cod. TAGASF32 H.3 - cod. TAGASF33 H.4 - cod. TAGASF34		cod. TAGCAP		 Cap basket cod. TAGCONT

meté
Biomedical

Via Boccaccio, 8
21010 Arsago Seprio (VA)
Tel./Fax: +39 0331.796417
www.dentalmete.it